




Silver Alginate Wound Dressing

Advanced Wound Care

We have harnessed the natural healing power of seaweed with this advanced wound care dressing. With high volume absorbtion, the dressing quickly absorbs exudates while forming a gel-like substrate covering over the wound. This substrate helps maintain a moist environment for ideal wound healing and keeps the dressing from adhering to the wound so it is painless to remove. Our silver particle additive creates an anti-bacterial shield that helps reduce the risk of infection.

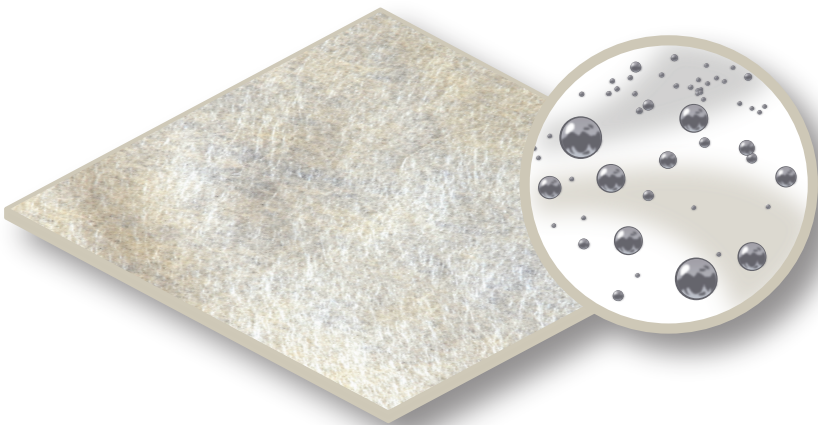


Usage

This type of dressing is best used when the risk of infection is high. Wounds appropriately treated with the dressing are moderately to heavily exuding partial to full thickness wounds including (but not limited to):

- Post-operative surgical wounds
- Graft and donor sites
- Trauma wounds
- Chemical burns
- Pressure ulcers
- Diabetic ulcers
- Leg ulcers

Ref Number	Description	QTY
9020104	2x2" (0.5cm x 0.5cm)	10/100
9020102	4x5" (10cm x 12.5cm)	10/100
9020103	4x8" (10cm x 20cm)	5/50
9020105	1x12" (0.25cm x 30cm)	5/50

For more info: 1.844.4HEALQU | www.healqu.com | info@healqu.com

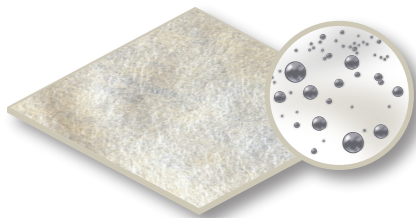

INTENDED USE:

Direct application to full and partial thickness, exuding wounds. This product provides excellent exudate absorption and maintains a moist, anti-microbial healing environment to promote wound healing.

INDICATIONS FOR USE:

Moderately to heavily exuding partial to full thickness wounds, including:

- Post-operative surgical wounds
- Graft and donor sites
- Trauma wounds
- Pressure ulcers
- Diabetic ulcers
- Leg ulcers



CONTRAINDICATIONS:

This product is not intended for use in the following circumstances:


- third-degree burns
- dry or lightly exuding wounds
- individuals with a known sensitivity to alginates or silver
- to control heavy bleeding
- surgical implantation

PRECAUTIONS


Systemic antimicrobial therapy should be considered when wound infection is evident. HEALQU™ Silver Alginate Dressing may be used, under medical supervision, in conjunction with systemic antibiotics. The dressing may adhere if used on dry or very lightly exuding wounds. If the dressing is not easily removed, moisten it with sterile saline solution prior to removal. The dressing performance may be impaired by excess use of petroleum-based ointments. For external use only.

DIRECTIONS FOR USE:


1. Clean the wound with sterile saline water and dry the surrounding skin gently.

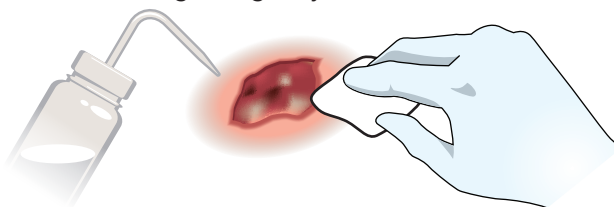

2. Remove dressing from packaging and apply to wound.

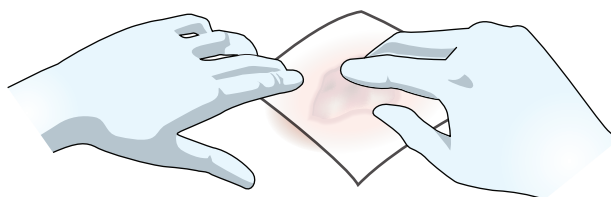

3. Cover and secure the HEALQU™ Silver Alginate Dressing with a secondary dressing.

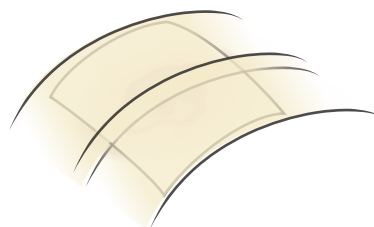

NOTES:

Dressing change frequency will depend on patient condition and the amount of exudates. Initially it may be necessary to change the dressing every 24 hours. When exudate quantity subsides, dressing change frequency may be reduced.

Observe the wound before removing dressing. If the wound appears dry, saturate the dressing with sterile saline water prior to removal.

